

Moxi Now at a glance

		Pro	800	700	600
Channels		20	20	16	10
Signature features	Glossary				
SpeechZone 2	In noisy situations, adjusts to an optimal microphone strategy to focus on speech from any direction	SpeechZone 2	SpeechZone		
Binaural spatial processing	Hearing instruments work together to determine the direction of speech	•	•		
SoundNav	Automatic program that identifies and smoothly transitions across environments, with a special focus on conversation types	7 environments	6 environments	5 environments	2 environments
Sound Conductor	Dynamically balances the appropriate features providing speech understanding and comfort while maintaining natural sound quality	•	•	•	•
Log It All	Captures time spent across all seven SoundNav environments to support technology recommendations	•	•	•	•
MyMusic	Specialized fitter-controlled settings to enhance music enjoyment	Automatic binaural	Automatic binaural	•	•
Binaural Phone	Streams conversation from phone ear to opposite ear, no accessory required	•	•	•	•
Automatic Adaptation Manager	Intelligently and gently transitions patients to the amplification required for long term optimal speech understanding from first fit acceptance	•	•	•	•
Features					
Adaptive directional	Directional strategy that tracks and suppresses multiple background noise sources	Multiband	Multiband	Multiband	Multiband
Frequency compression	Compresses and shifts inaudible high frequency sounds into a more audible lower frequency region	•	•	•	•
AntiShock	Eliminates discomfort from sudden impulse noises	•	•	•	•

In all technology levels:

- Pinna Effect
- Natural Sound Balance
- Data logging
- Feedback manager
- Wind control
- Tinnitus masker
- DuoLink
- Easy-t
- Plasma coating
- WDRC processing with NAL-NL1, NAL-NL2 and DSL v5
- Linear processing with DSL v5


Accessories give you another tool to solve patient challenges. They are available at all technology levels.

At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to provide hearing solutions that improve lives in meaningful ways. Because hearing matters.

© 2016 Unitron. All rights reserved.

unitron.com

16-027 028-6229-02

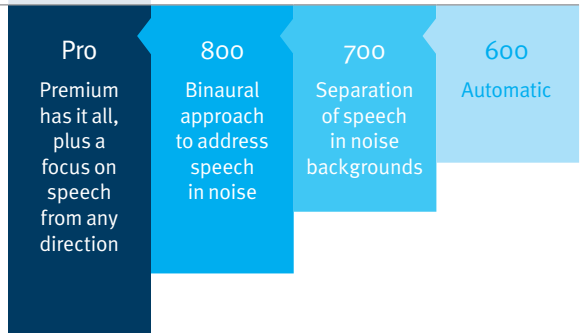

Smallest size doesn't mean compromise

It may be the smallest, but Moxi Now is packed with amazing technologies that will make patients really happy.

sonova
HEAR THE WORLD